

Overview: Omnilux Contour

Your best skin starts here.

Our clinically-proven Contour line of LED products utilize red (633nm) and near-infrared (830nm) light to boost collagen production and soothe skin in just 10 minutes.





Contour Face

Dermatologist recommended red and near-infrared light therapy for clinically proven results for anti-aging and skin rejuvenation.



Contour Neck & Décolleté

Award-winning red and near-infrared light therapy device with clinically-proven results for anti-aging and skin rejuvenation.



Contour Glove

Quick and effective red and near-infrared light therapy with clinically-proven results for skin rejuvenation and joint pain.

Contour Face



Clinically proven benefits:


1. Reduces fine lines and wrinkles
2. Reduces the appearance of pigmentation and redness
3. Promotes healthier, younger-looking skin
4. All natural. No pain, side effects, or downtime



Protect your most recognizable feature


Omnilux Contour Face is the most powerful flexible LED device available. Featuring 132 medical-grade LEDs, and optimized wavelengths of red (633nm) and near-infrared (830nm) light, Contour Face is designed and manufactured using Omnilux's proprietary, gold standard medical device technology - providing you derm-level results right from home.

- ✓ Stimulates collagen production for firmer, more youthful-looking skin
- ✓ Refines skin texture and tone for an effortless glow
- ✓ A non-invasive and natural option
- ✓ Maintains and improves in-clinic treatments

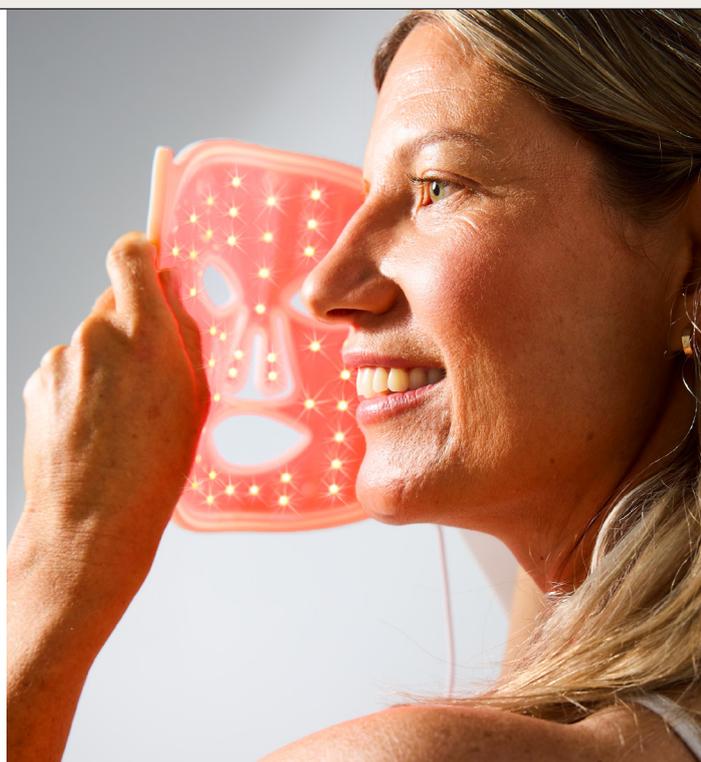

Clinically proven results:

After 3 treatments over 4 weeks, users agreed:

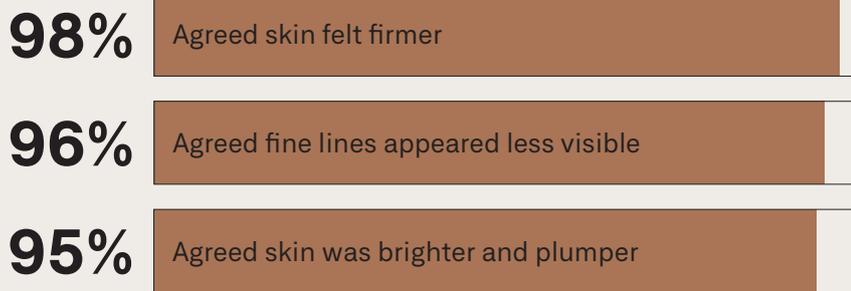

Technical specifications:

<p>LED TYPE</p> <p>Red and near-infrared</p>	<p>TARGETED USE</p> <p>Face</p>	<p>INCLUDES</p> <ul style="list-style-type: none"> • Flexible silicone Face device • Rechargeable battery controller • Power supply and plug adaptors • USB-C line to USB-A connector • 2x adjustable straps
<p>WAVELENGTHS</p> <p>633nm (red) and 830nm (near-infrared)</p>	<p>CERTIFICATIONS</p> <p>FDA, CE, TGA</p>	
<p>NUMBER OF LEDs</p> <p>132 LEDs (2 per bulb)</p> <p>66 LED bulbs</p>	<p>WARRANTY</p> <p>2 years: device</p> <p>2 years: controller</p>	



How to use:

We recommend three to five 10-minute treatments per week for four to six weeks, followed by maintenance as desired.

- Charge your controller & then connect it to the device
- Prepare your skin by cleansing and patting dry
- Secure the device to your head using the velcro straps
- Start the 10-minute treatment
- Remove the mask and apply moisturizer to the skin
- Gently wipe down the mask

Contour Neck and Décolleté



Clinically proven benefits:

1. Reduces fine lines and wrinkles on the neck and chest
2. Reduces redness, inflammation, and the appearance of sun damage
3. Promotes healthier, younger-looking skin
4. All natural. No pain, side effects, or downtime



Treat beyond the face

Too often people neglect their neck and décolleté. This area is much thinner than facial skin and requires just as much attention. Our Contour Neck & Décolleté device offers a dynamic design that allows you to optimally target these overlooked treatment areas. Perfect partner to the Contour Face mask makes for the ultimate skincare experience.



Real users, real results



Technical specifications:

LED TYPE

Red and near-infrared

WAVELENGTHS

633nm (red) and 830nm (near-infrared)

NUMBER OF LEDs

108 LEDs (2 per LED bulb)

54 LED bulbs

TARGETED USE

Neck and décolleté (upper chest)

CERTIFICATIONS

FDA, CE, TGA


WARRANTY

2 years: device

2 years: controller

INCLUDES

- Flexible silicone Neck & Décolleté device
- Rechargeable controller
- Power supply & plug adaptors
- USB C line to USB A connector
- 2x adjustable straps

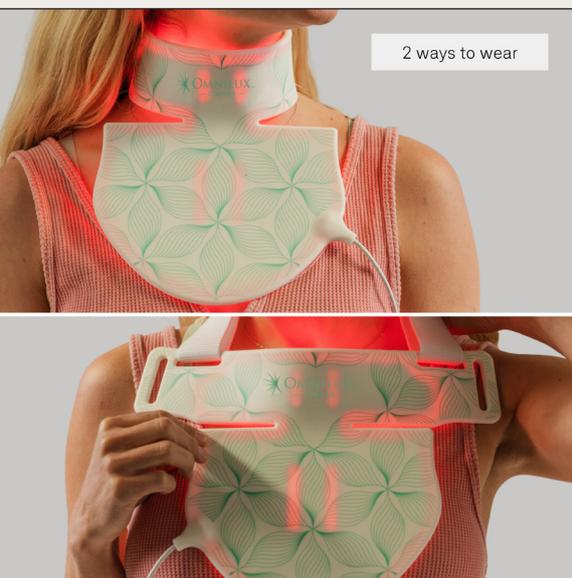

How to use:

We recommend three to five 10-minute treatments per week for four to six weeks, followed by maintenance as desired.

- Charge your controller & then connect it to the device
- Prepare your skin by cleansing and patting dry
- Secure the device around your neck using the velcro straps
- Start the 10-minute treatment
- Remove the device and apply skincare to the area
- Gently wipe down the device

Contour Glove



Why Contour Glove?

Omnilux's therapeutic wavelengths of light penetrate deep into the skin, promoting cellular regeneration and collagen production to foster both external changes (skin rejuvenation) and internal healing (reduced inflammation).

Equipped with 144 LEDs and manufactured using Omnilux's medical-grade LEDs, the Contour Glove provides a safe and effective treatment built to protect and correct an often overworked and under treated area.

How to use:

We recommend three to five 10-minute treatments per week for four to six weeks, followed by maintenance as desired.

- Charge your device
- Secure the device to your hand and wrist using the velcro strap
- Start the 10-minute treatment
- Remove the device and apply moisturizer to skin
- Gently wipe down the device



Clinically proven benefits:

1. Improves the appearance of sunspots, pigmentation, wrinkles, dry cracking skin, and uneven skin tone
2. Reduces inflammation, pain, and irritation associated with aging joints*
3. Promotes healthier, younger-looking skin
4. All natural. No pain, side effects, or downtime

*Omnilux Contour Glove is regulated by the FDA for the treatment of pain. Outside of the US, Omnilux Contour Glove is pending approval for pain management.



Real users, real results



Technical specifications:

LED TYPE

Red and near-infrared

WAVELENGTHS

633nm (red) and 830nm (near-infrared)

NUMBER OF LEDS

144 LEDs (2 per LED bulb)
72 LED bulbs

TARGETED USE

Hand and wrist

CERTIFICATIONS


FDA, CE, TGA

WARRANTY

2 years: device

INCLUDES

- Flexible silicone Glove device
- Power supply & plug adaptors
- USB C line to outlet connector
- Adjustable strap
- Carry bag and built-in battery



Ask your provider how to incorporate Omnilux
Contour Line into your routine, or visit

www.omniluxLED.com

+1-707-309-7024

support@omniluxLED.com

