

Omnilux Men



Engineered for men's skin

Omnilux Men is the only LED mask specifically engineered for men's specific skin needs.

From razor burn & dryness to fine lines, dark circles and under-eye bags, Omnilux's advanced LED light therapy technology provides visible, dermatologist-level results proven for men.



Level up with 1072nm

Equipped with red and near-infrared optimized LEDs, Omnilux Men incorporates a powerful 1072nm near-infrared wavelength of light. This longer beam of light targets fibroblasts deep within the skin, delivering the most effective treatment for men.

1. Perfectly activates cellular structures in thicker skin
2. Reduces fine lines, wrinkles, redness, irritation, dark circles, & under-eye bags
3. Promotes a smoother, more even-toned complexion



How to use:

We recommend three to five 10-minute treatments per week for four to six weeks, followed by maintenance as desired.

- Charge your controller & then connect it to the device
- Prepare your skin by cleansing and patting dry
- Secure the device to your head using the velcro straps
- Start the 10-minute treatment
- Remove the mask and apply moisturizer to the skin
- Gently wipe down the mask

Real users, real results



Technical specifications:

LED TYPE

Red and near-infrared

WAVELENGTHS

633nm (red)

830nm (near-infrared)

1072nm (near-infrared)

NUMBER OF LEDS

132 LEDs (2 per LED bulb)

66 LED bulbs

CERTIFICATIONS

FDA, CE, TGA

WARRANTY

2 years: device

2 years: controller


INCLUDES

- Flexible silicone Omnilux Men device
- Rechargeable battery controller
- Power supply and plug adaptors
- USB-C line to USB-A connector
- 2x adjustable straps
- Carry bag

In clinical studies, Users agreed:

After 3-5 treatments per week for six weeks

- ✓ Reduction in UV spots (sun damage)
- ✓ Fine lines appeared less visible
- ✓ Skin tone was more even

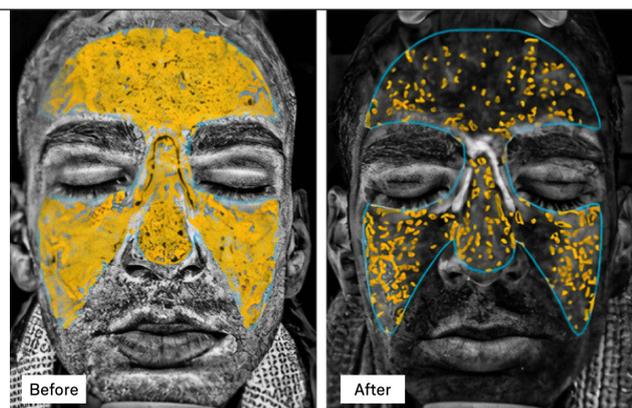

Ask your provider how to incorporate Omnilux Men into your routine, or visit
www.omniluxLED.com

+1-707-309-7024

support@omniluxLED.com

